

Cuticell® Contact®

Silicone Wound Contact Layer



Highly Flexible

Atraumatic

**Single-sided
adherence**

Cuticell® Contact – impressively **wound calming**

Cuticell® Contact is a silicone wound contact layer that has a perforated transparent structure and is highly conformable. This makes **Cuticell® Contact** an effective primary wound dressing for advanced wound management.

Wound-friendly right from the start

Once applied, **Cuticell® Contact** can remain in place for up to 14 days, protecting the wound. Exudate can easily pass through the perforations and there is minimal risk of tissue in-growth. As the adhesive is only on the wound contact side, secondary dressings can be easily changed without causing trauma to the wound.

- Undisturbed healing
- Pain-free, atraumatic dressing changes
- Improved patient compliance
- Cost-effective



Impressive technology – exceptional **security**

Crystal-clear transparency

Cuticell® Contact's high transparency simplifies application.

- Precise application
- Simple wound inspection without removal of the silicone dressing

Optimal perforated structure

With up to ten perforations per cm² **Cuticell® Contact** ensures even highly viscous exudate can easily pass through.

- For wounds with low to high exudate levels
- Helps prevent maceration

Unique elasticity

The combination of a thin PU film and a one-sided silicone coating make **Cuticell® Contact** highly flexible and conformable.

- Conforms to difficult to dress areas
- Stays in place during secondary dressing change



Siltec Technology

Cuticell® Contact benefits from **Siltec Technology**. The silicone-coated wound contact side only adheres to dry healthy skin, therefore offering very gentle adhesion.


- Optimal protection of skin and wound
- Secure adhesion without sticking to wound bed
- Tissue-friendly and pain-free removal
- Gentle removal without residue
- Hypoallergenic

Easy to apply – unique in performance




Step 1

Remove the overlapping protection film from the adhesive face of the dressing.

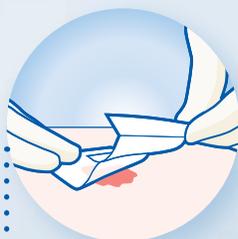

Step 2

With the adhesive face down, apply **Cuticell® Contact** making sure to overlap wound margin by 2 cm (3/4 in).

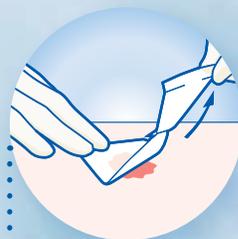

Step 3

Now remove the second protective film from the dressing surface and smooth softly with your hand.

Cuticell® Contact may be used in conjunction with a secondary dressing based on drainage and secondary fixation.

- Absorbent Dressing
- Cutisorb® Aginate

For fixation purposes the following products are recommended:

- Elastomull®
- Cover-Roll® / Cover-Roll® Stretch
- Fixomull® transparent

Cuticell® Contact – indications

Silicone wound contact layers are suitable for use on the following wound types:

- Cuts
 - Abrasions
 - Burns
 - Chronic wounds (also under compression)
 - Lacerations
 - Blisters
 - Skin lesions
 - Also appropriate for skin protection (e.g. after radiation treatment)
-
- Can be used in combination with antiseptics
 - Suitable for very sensitive and fragile skin

All benefits at a glance

Cuticell® Contact impresses users and patients with a multitude of benefits:

- Longevity and cost effective
- Atraumatic dressing changes
- No sticking to gloves
- No exudate congestion
- Minimised risk of maceration
- Easy wound inspection
- Flexible and conformable
- Easy application



Cuticell® Contact®

Cuticell® Contact is indicated for the following wound conditions:

Wound depth  Shallow

Exudate level  Slight to moderate

Wound phase  Fibrin coating

 Granulation

 Epithelialisation

 Moderate to strong

 Strong to excessive



Ordering Information:

REF No.	Size	Dressings/Box	Box/Shipper	HCPCS*
7268000	2in x 3in (5 cm x 7.5 cm)	5	30	A6206
7268001	3in x 4in (7.5cm x 10 cm)	5	20	A6206
7268002	4in x 7in (10 cm x 18 cm)	5	20	A6207
7268003	6in x 10in (15 cm x 25 cm)	5	10	A6208

Please refer to the product label and/or package insert for full instructions on the safe use of these products.

* This document includes codes that might be used to bill for BSN medical products and related services. It is the provider's sole responsibility to determine and submit appropriate codes, charges, and modifiers for services rendered. BSN medical Inc. cannot guarantee coverage or reimbursement.

BSN medical Inc. 5825 Carnegie Blvd. Charlotte, NC 28209-4633
Tel. (+1) 704 554 9933 Fax (+1) 704 358 4558 www.bsnmedical.com
To order toll-free: BSN medical (+1) 800 552 1157

