

Cutimed®

HydroControl

Cutimed®
HydroControl

One dressing
for both dry
and moderately
exuding
wounds




Closing wounds. Together.

 **essity**

Cutimed® HydroControl: A new way to manage wounds

INNOVATIVE 4 - PRONGED APPROACH

While many dressings focus on just one element of wound management, **Cutimed® HydroControl** is based on an innovative product concept uniquely combining four key areas to promote successful wound healing.

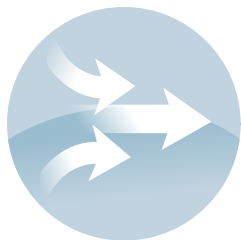


HIGH WOUND STIMULATION POTENTIAL

Stimulates and supports autolytic debridement by an osmotic effect.




EXCELLENT EXUDATE MANAGEMENT

Promotes high therapeutic safety as fluid is absorbed vertically and retained within the dressing.

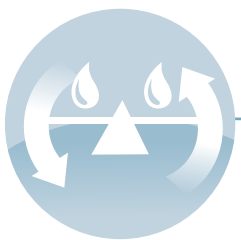

OPTIMUM MOISTURE BALANCE

The combination of Flexipores, extra-thin PU film and a hydropolymer gel matrix provides an ideal environment to support moist wound healing.

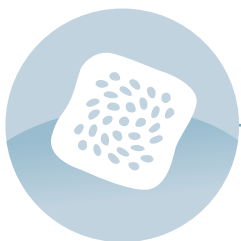


INNOVATIVE DESIGN

Thin, cuttable and easy to apply, the dressing combines patient comfort with high product performance.

• HYDROPOLYMER GEL MATRIX

• EXTRA-THIN PU FILM

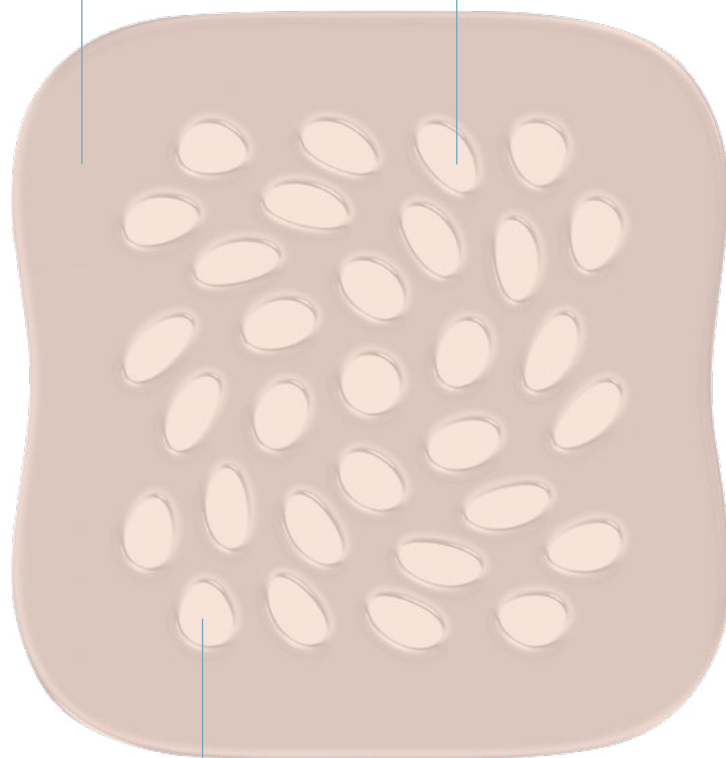


• FLEXIPORES

Wound stimulation – The key to successfully closing wounds

CREATING IDEAL HEALING CONDITIONS WITH A PHYSICAL MODE OF ACTION

Cutimed® HydroControl works by the physical principle of osmosis. As the dressing absorbs exudate, the osmolarity of the wound bed increases. This stimulates the perfusion of fresh fluid into the wound bed from surrounding tissues, supporting autolytic debridement. This makes **Cutimed® HydroControl** suitable for many different wound types and particularly for the use on stagnant wounds.

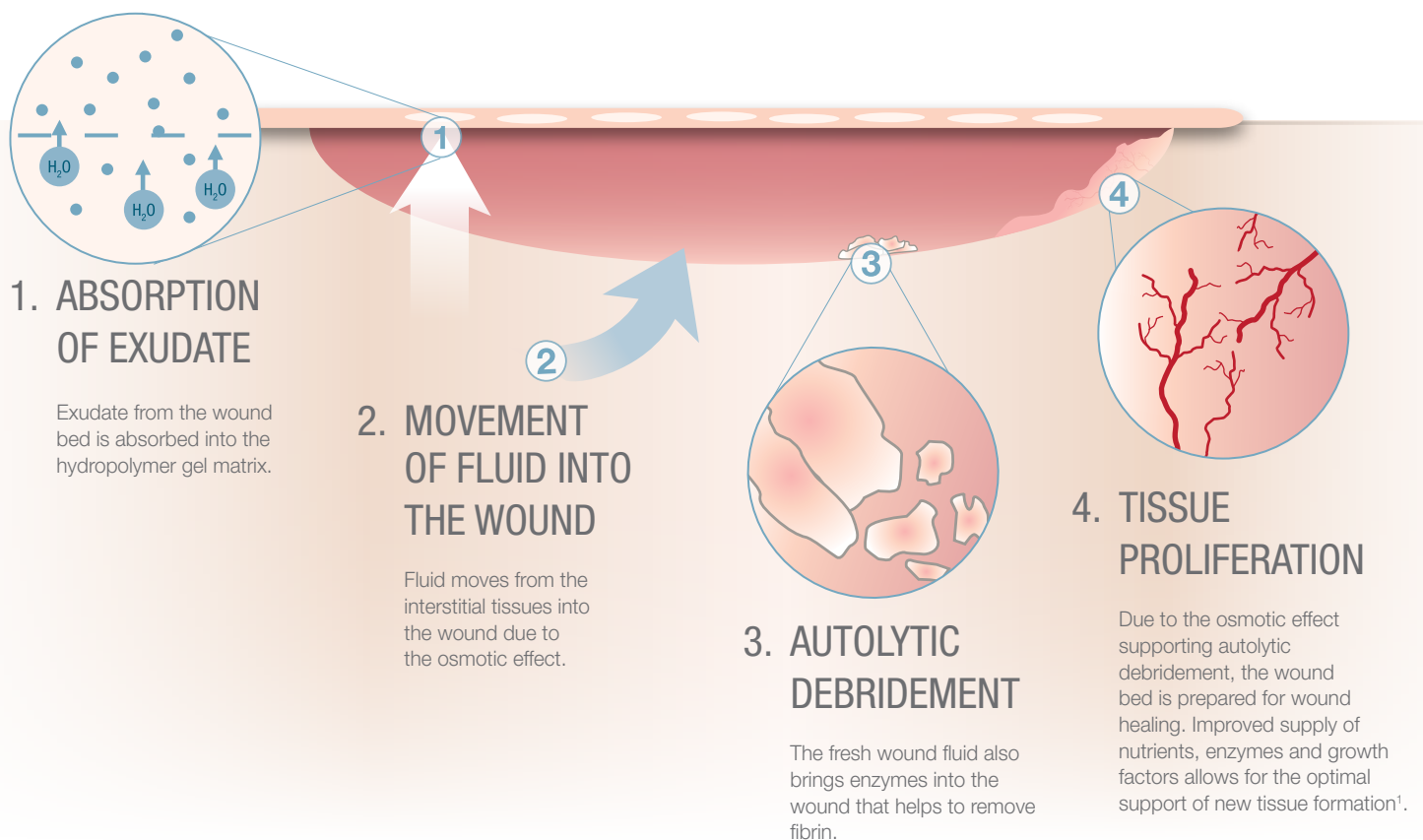


Ideal Moisture Management

IN EVERY WAY SUPERIOR

No other dressing has the same excellent combination of fluid management and moisture donation as **Cutimed® HydroControl**. In comparison to wound dressings for light to moderately exuding wounds, **Cutimed® HydroControl** demonstrates superior fluid handling due to its innovative composition.



Excellent results*

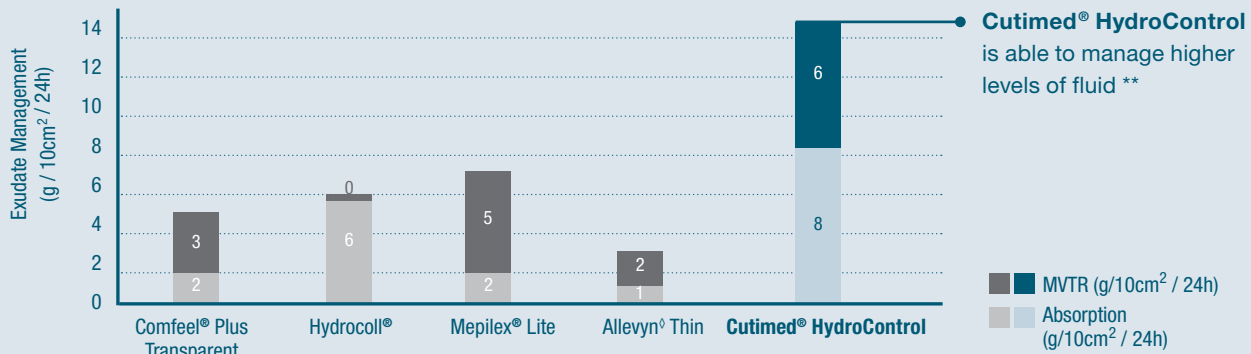

Cutimed® is a registered trademark of BSN medical GmbH. Comfeel® Plus Transparent is the trademark of Coloplast A/S, Mepilex® Lite is the trademark of Mölnlycke Health Care AB, ALLEVYN® Thin is the trademark of T.J. Smith & Nephew Ltd. and Hydrocoll® Border is the trademark of Paul Hartmann AG. BSN medical GmbH is not associated with any of these companies regarding the trademarks and respective products.

Laboratory testing results are intended to illustrate a product's performance under controlled conditions. Actual use results may vary.

Indications

Cutimed® HydroControl is indicated for dry to moderately exuding, shallow, chronic and secondary healing wounds such as:

- Venous or arterial leg ulcers
- Pressure ulcers
- Diabetic foot ulcers

Cutimed® HydroControl can be used under compression.

Cutimed® HydroControl

Size	Ref-No.	Dressings / Box	Suggested HCPCD Code
1¾ x 1¾ in	7323100	10	A6234
3 x 3 in	7323101	10	A6234
4 x 4 in	7323102	10	A6234
6 x 6 in	7323103	5	A6235

1 Weir D. How To... Top tips for wound dressing selection. Wounds International 2012, 3(4): 18-22

* Comparison: Cutimed® HydroControl compared to other products for low- to moderately exuding wounds.

** Source: BSN medical data in file report. JB4 40022569 - fluid handling; JB4 40022570 -vertical absorption

REIMBURSEMENT DISCLAIMER

This document includes codes that might be used to bill for BSN medical Inc. products and related services. Each provider will have to verify the appropriate codes for each patient. It is the provider's sole responsibility to determine and submit appropriate codes, charges, and modifier for services rendered. Providers should contact insurers to verify correct coding procedures prior to submitting claims related to any product or device. BSN medical Inc. cannot guarantee coverage or reimbursement with the codes listed in this billing guide. In all cases, providers will need to follow local payer policies for billing and reimbursement.

Please refer to the product label and/or package insert for full instructions on the safe use of these products.